

Chemical Terrorism Packaging and Shipping Flow Chart


STATE LABORATORY INSTITUTE

305 South Street

Jamaica Plain, Ma, 02130

1

- Use evidence tape over top of primary container.
- Initial half on the tape and half on the container.
- Treat blanks in the same manner

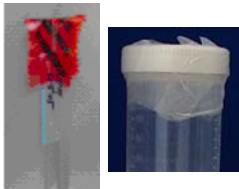

2

Pack primary containers in a way that prevents contact between tubes.



3

- Place wrapped primary containers inside leak-proof secondary container.
- Add enough absorbent to absorb the entire contents.



4

- Place evidence tape over the closure of the secondary container.
- Place initial half on the tape and half on the packaging.
- Place Biohazard Label on Secondary Package
- Place Specimen Submission forms in pouch of secondary container



5

Place secondary package in outer container.



6

- Use a Styrofoam insulated shipper.
- Add cushioning and **cold packs for blood**, or **dry ice for urine** to the bottom



7

- Add more cushioning material.
- Place **cold packs** or **dry ice** on top.
- Place shipping manifest inside outer box, wrapped in a zip top bag.
- Chain of custody should not be included



8

Label boxes accordingly



Labels for Blood Shipment:

To (shipper) From (consignee)
Diagnostic Specimens UN3373
Orientation arrows

Labels for Urine Shipment:

To (shipper) From (consignee)
Diagnostic Specimens UN3373
Orientation arrows Class 9
UN 1845